

What is so special about Aesculap clipper products?

Aesculap Clippers are:

- very robust and reliable
- can manage even difficult coats
- high-quality components and processing
- Made in Germany

...and the cutter heads and shearing blades have the following advantages:

- the best shearing properties
- long-lasting cutting performance
- can be reground up to 7 times
- the development and production of the machines and cutter heads at Aesculap Schermaschinen GmbH in Suhl
- expertise in developing clippers, cutter heads and shearing blades

Aesculap clippers – quality from Germany – since 1912

Got questions about our clippers?
Our product experts will be pleased to answer all of your questions about the types of machines, shearing blades and the professional use of the equipment.

Germany hotline:
+49 8086 933-552

Hotline in UK:
+44 (0) 1572 722558

Email:
schermaschinen@kerbl.com

E-mail :
enquiries@kerbl.co.uk

Albert Kerbl GmbH
Felizenzell 9
84428 Buchbach
Germany

Kerbl UK Limited
8 Lands End Way
Oakham, Rutland
LE15 6RF UK


www.kerbl.com

www.kerbl.co.uk

Overview of selected Aesculap A5 type SnapOn blade sets for VET requirements

- very long edge-retention thanks to high-quality carbon steel
- quality manufacture at CNC machining centres by Aesculap Schermaschinen GmbH
- all parts are milled or profile ground, giving them a high level of stability
- high-end surface treatment ensures high level of corrosion protection

You can find more SnapOn shaving heads on our website.

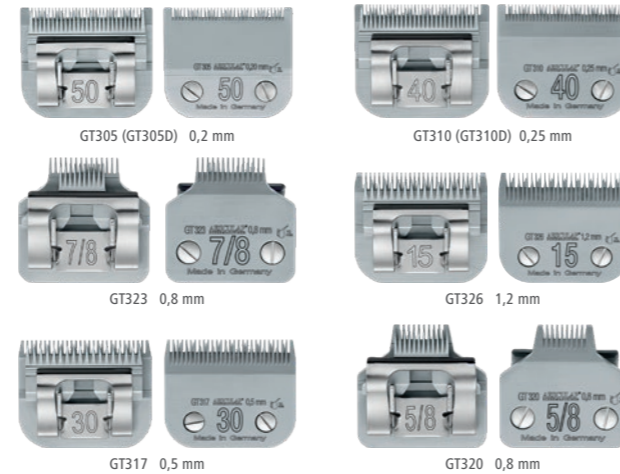

NEW Blade Sets with DLC Coating

The Aesculap "Diamond-like Carbon" (DLC) coating perfectly combines the property profiles of diamond- and Teflon-based surface enhancement.

This special coating guarantees outstanding product properties:

- high wear protection and longer service life
- excellent gliding thanks to minimised frictional resistance
- extended edge retention
- increased corrosion protection
- reduced care requirement thanks to dirt-repellent surface
- less heat development

GT330 1,5 mm GT345 3,2 mm GT360 6,3 mm GT364 9,5 mm

Changing the cutter head

Removing the cutter head:

1. loosen the cutter head lock by pressing the push button
2. press against the cutter head from the front to tip it backwards
3. pull the cutter head out of the blade hinge



Attaching the cutter head:

4. place the cutter head on the blade hinge
5. switch on the machine
6. press the cutter head forwards until it clicks



Important: the machine must be turned on so that the rocker arm can snap into the cutter head. If it does not snap in, you will hear a "rattling" sound and the cutter head will not cut.



Tip: If the blade hinge for the cutter head is accidentally tipped, you can loosen it again using a screwdriver.



Cleaning and caring for the cutter head:

GT604 Aesculap Oil
GTA103 Aesculap Blade Cool 2.0



We recommend oiling the cutter head before each shearing session. You can both reduce heat build-up and clean at the same time with the Blade Cool. Use this at least every 10 minutes. Oiling during shearing also improves the cut and extends the sharpness of the cutter head.



AESCULAP®
CLIPPERS

Quality from Germany – since 1912



Durati Vet

Professional battery-powered clippers for veterinarians

Aesculap Durati Vet

cordless clippers with Removable Battery for SnapOn Cutter Heads

The Durati cordless clippers are made with the proven Aesculap quality. These universal clippers with their standard SnapOn cutter head system (compatible with manufacturers like Heiniger, Oster, Wahl and Andis) boast an unbeatable price/performance ratio. With a balanced weight of just 380 g including battery and cutter head, and its lithium-ion battery that enables up to 4 hours of cutting, the Durati provides low-fatigue shearing. A second battery is optionally available and allows continuous operation. Effortless shearing with no annoying cables!

- best value for money
- Aesculap Performance
- for veterinary care of dogs and small animals

Features:

- Battery run time: up to 4 hours
- Battery charge time: 4 hours
- Battery capacity: 2850 mAh
- new Li-ion battery technology (fast charge and durable)
- lightweight – just 380 g including cutter head and battery
- very quiet – noise level only 63 dB

An optional and extensive SnapOn cutter head range means you can shear all kinds of coat types and up to 9.5 mm.

Set GT454-NB includes:

- 1 x cordless clippers Durati VET
- 1 x battery
- 1 x charging station with universal plug for various countries
- 1 x cutter head GT310 #40 (0.25 mm)
- 1 x Aesculap oil (GT604)
- 1 x instruction manual



Aesculap SnapOn Blade Sets

Precision from the profi, for the profi!

Blade set from Aesculap Schermaschinen are precision parts which aim at perfection.

The blade set determines the performance of every clipper. The SnapOn blade set from the metal forge in Suhl/Thuringia are genuine "Made-in-Germany" products.

Exceptionally high-quality cutting tools are manufactured from top-quality carbon blanks in several production stages, in compliance with medical technology directives. Except for the hardening process, the parts never leave the production facility.

This in-house production can only be to your advantage:

- **unmatched long cutting durability**, thanks to first-grade carbon steel
- **consistently high cutting capacity** of each individual blade set, thanks to state-of-the-art CNC grinding and milling equipment
- **exceptional coat gliding behaviour**, thanks to high-end surface finish
- **uniform and clean cutting pattern**, thanks to the specially developed tooth geometry
- **extremely regrindable, often up to 7 times**, thanks to the combination of material hardness, quality and manufacturing expertise

precise and consistent plate pressure for controlled heat generation



Precision hollow-ground for the best cutting results

specially developed tooth geometry for precise immersion in the coat



exceptional coat gliding behaviour due to high-end surface finish

first-class raw materials for long-lasting cutting, and up to 7 regrinds



Which cutter head should I use for which coat?

The cutter head advisor for SnapOn cutter heads

It is important to use the right cutter head for each coat type and for the various parts of an animal. There is a wide range of Aesculap SnapOn cutter heads available.

These figures are indicative and based on experience; coat structures can vary by individual animal. We recommend analysing the coat structure before each use of the clippers.

Cutter head	Designation	Shear head size	Shearing height	Shear width	Use / application for animal breeds	Coat structure	Body part	VET application
GT305 GT305D	Cutter Head SnapOn 0.2 mm	50	0,2 mm	40 mm	mixed breeds, cats, rabbits	all coat structures	preparation of surgical sites	venupuncture, preparation for surgery
GT310 GT310D	Cutter Head SnapOn 0.25 mm	40	0,25 mm	40 mm	mixed breeds, cats, rabbits	all coat structures	preparation of surgical sites	biopsy areas, preparation for surgery
GT317	Cutter Head SnapOn 0.5 mm	30	0,5 mm	40 mm	mixed breeds, cats, rabbits	matted coat	Body	Ultrasound preparation
GT320	Cutter Head SnapOn 0.8 mm	5/8	0,8 mm	20 mm	all breeds		Paw shearing	
GT323	Cutter Head SnapOn 0.8 mm	7/8	0,8 mm	20 mm	all breeds		Paw shearing	
GT326	Cutter Head SnapOn 1.2 mm	15	1,2 mm	40 mm	mixed breeds, all cut breeds	matted coat	Body	Ultrasound preparation