#11602

AniShock Comfy Pro

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(EN) Operating instructions for the AniShock Comfy Pro (ref. no. 11602)

safety instructions
Please read the following notes carefully before using the device.

CAUTION!
High voltage! Danger of personal injury!

- Do not touch the electrodes. A voltage pulse may still be generatedeven when the device is switched off (residual charge).
- Never use the device on persons.
- Keep the device out of the reach of children.
- Remove the batteries and fully discharge the device before starting the cleaning process.
- Do not operatethe device whilst it is moist.
- Do not operate the device in liquids.
- Do not use the device for striking.


CAUTION!
Danger of personal injury and damage to property due to improper use!

- Do not expose the device to excessive temperatures, direct sunlight or fire!
- Only charge the battery using the power pack supplied with the original device!
- Do not leave the battery charging unattended.
- Stop charging immediately if you notice the battery becoming excessively hot
 .
- Do not use the device if the casing is broken!

Intended use
The device is intended as a herding aid for fully-grown cattle and pigs. Any other use, in particular applications on persons, is considered improper. The animal welfare regulations applicable in the respective country must be observed.

Function
The device generates a high voltage and delivers an electric shock when the animals come into contact with the electrodes.

Construction
1 electrode
2 Shaft
3 LED luminaire
4 LED buttons
5 Impulse button
6 Wrist strap

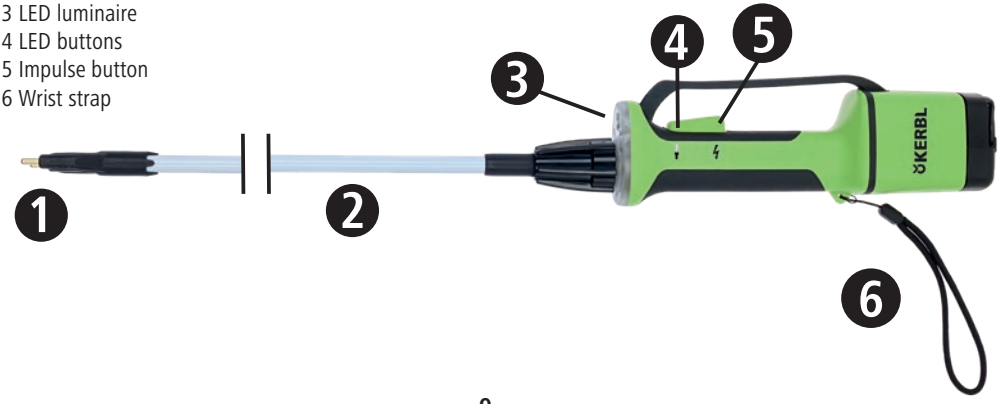7 On/off switch with battery indicator
8 Charging indicator
9 Charging cable connection

Charging the battery
During operation, the battery indicator on the on/off switch (7) shows the battery level:
- Green: Sufficient battery charge
- Red: Battery level below 20% – charge as soon as possible
- Flashing red: Battery critically low – charge immediately

Only use the supplied charging cable and power pack for charging.

1. Switch off the device (7).
2. Remove the protective cap from the charging port (9) on the battery pack. Ensure that the contact surfaces are clean and dry.
3. Connect the charging cable with the magnetic connector to the charging port.
4. Connect the charging cable to a USB socket or to the mains (230 V) via the supplied power pack.
5. Once charging is complete, disconnect the charging cable and power pack, and replace the protective cap on the charging port.

During charging, the charge indicator (8) shows the charging status:
- Red: Battery is charging
- Green: Battery fully charged

If the charging indicator is not lit, check the magnetic contact and the power supply.

The charging time is max. 5.5 hours.

Operation

CAUTION! After pressing the impulse button, the device may remain electrically charged for up to 20 seconds. Do not touch the electrodes and avoid accidental contact with people or animals.

1. Switch on the device (7).
2. Place the electrodes (1) on the animals' hind leg muscles and press the impulse button (5). Whilst the button is pressed – but for no longer than one second – an electrical impulse and an audible signal are emitted.
3. Remove the device from the animal and switch it off.

Tip – for the gentlest possible use with an audible signal:
- Press the impulse button short and release it without touching the animal with the electrodes.
- If the animal does not react to the audible signal, attach the electrodes to the animal within 20 seconds. The stored electrical charge is then released as a single, short impulse.

Operating the LED luminaire
1. Switch on the device (7).
2. Press the LED button (4) to switch the LED luminaire on or off.

Replacing the battery pack
Use only original battery packs.
1. Switch off the device.
2. Pull the battery pack out of the device using the tab.
3. Insert the new battery pack fully into the device and ensure it clicks securely into place.

cleaning

CAUTION! High voltage!
Switch off the device and discharge any remaining charge before cleaning. Allow the device to dry completely before using it for further use. The device can be cleaned using standard cleaning agents and water.

Disposal

The device must not be disposed of with household waste. It must be taken to a suitable collection or return point for waste electrical and electronic equipment. Before disposal, remove the battery pack from the device and dispose of it separately.

The battery must not be disposed of with household waste. Used batteries must be collected separately and taken to a suitable collection or return point.

technical data

Output voltage:	max. 10 kV / average 6 kV
Output current:	≤ 1 mA
Max. pulse duration:	1 s
Battery:	Li-ion, 7.2 V / 2,600 mAh / 18.72 Wh
Mains adapter input:	100–240 V AC, 50/60 Hz
Charging output:	5 V DC / 2 A
Charging time:	approx. 5.5 hours
Protection type:	IP 45

spare parts

Ref. no.	designation
11607	Shaft, 71 cm
11608	shaft, 58 cm
11611	charging cable and power pack
11613	battery pack

warranty
This product is covered by a 24-month warranty from the date of purchase. For further information and warranty terms and conditions, visit kerbl.com/guarantee.

EU Conformity
Albert Kerbl GmbH hereby declares that the product described in this manual complies with the relevant EU harmonisation legislation.

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AB Uygunluğu
Albert Kerbl GmbH, bu kılavuzda açıklanan ürünün ilgili AB uyum mevzuatına şuna uygundur dediğini beyan eder.
Kerbl UK Ltd